**AQUATIS**

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### Installing the batteries

- Open the battery compartment
- Install the batteries in the correct direction.
- Close the battery compartment


### Install the APP

- Activate Bluetooth on your Smartphone
- Your Smartphone and tablet must be connected to the network (WiFi, 4G or 3G)
- Use your device to scan the QR code or search for "OKOK International" in App Store or Google Play to install the application.

### Indications for use

Register and create a new user, if you decide to register with an email, follow the instructions below. Select the option "Register with email", enter your email and press the "Next" button

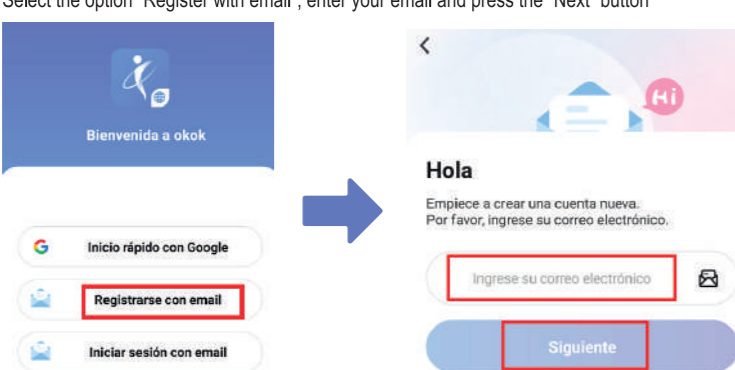
You will receive a code in the inbox of the email you entered above (if you do not receive the email please check your Spam folder). Once you get the code, enter it and press the "Next" button


button, then it will ask you to reset your password, enter one you can remember, accept the terms and conditions and press the "Register" button.


Fill in your personal information and click the "Next" button.

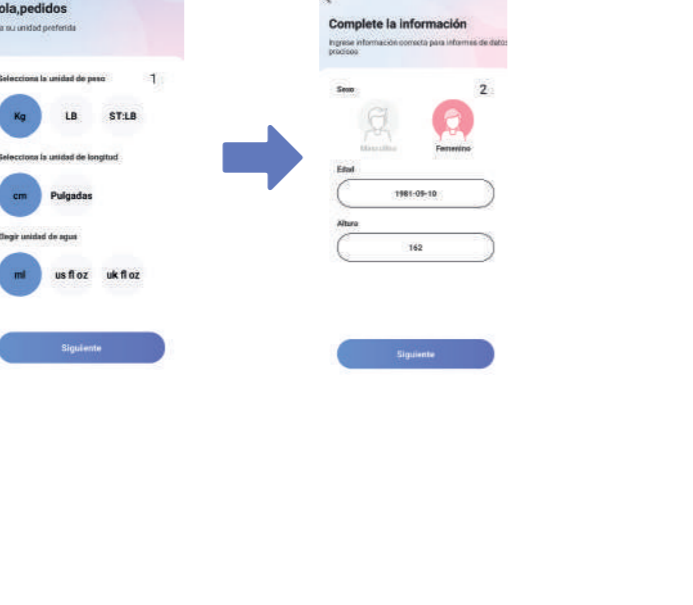Link the scale by pressing the "Add device" button

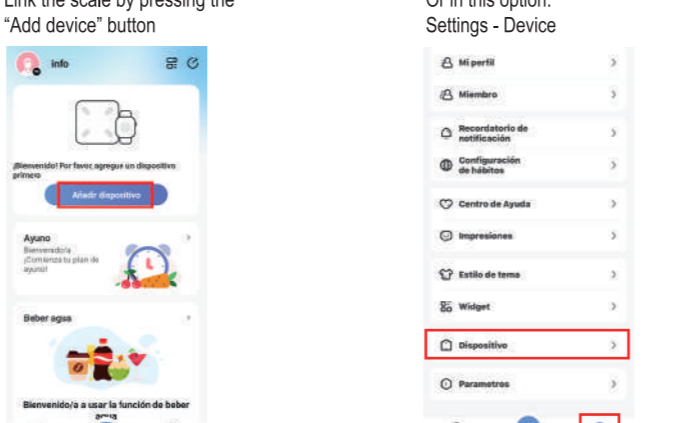
Or in this option: Settings - Device

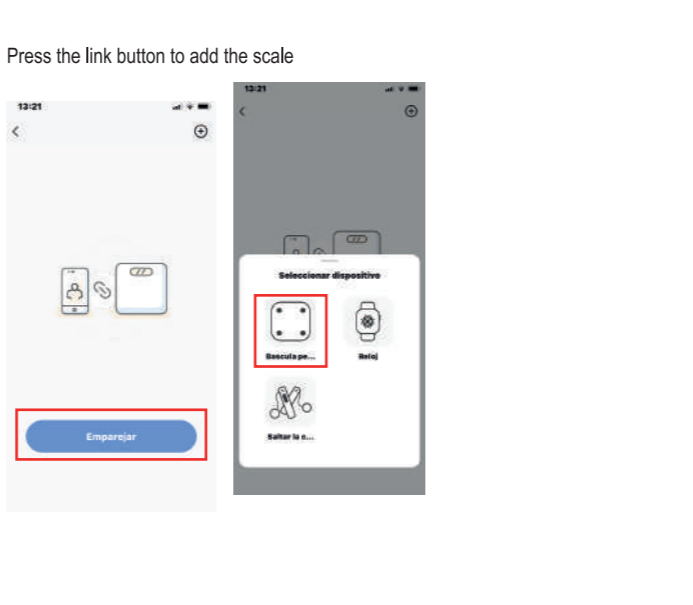
Press the link button to add the scale


Active all Bluetooth permissions and allow your mobile device to search for devices

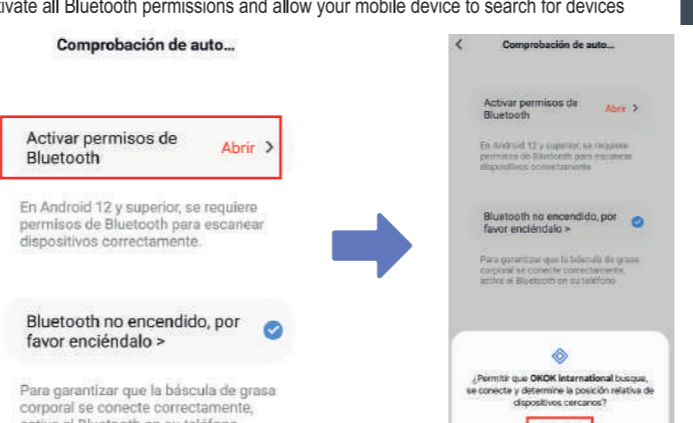
Go up the scale so the app can detect you and record your first weight. This will take a few seconds.

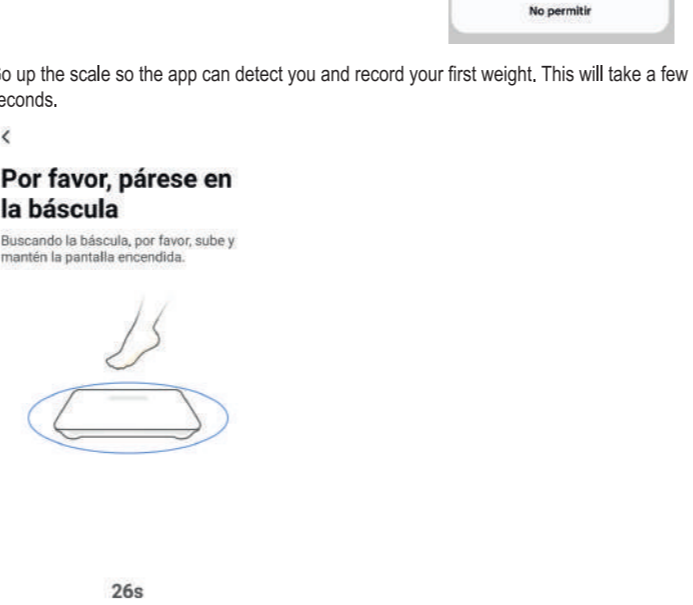If you want to weigh yourself again, press the button , go up to the scale and press the "Start Weighing" button

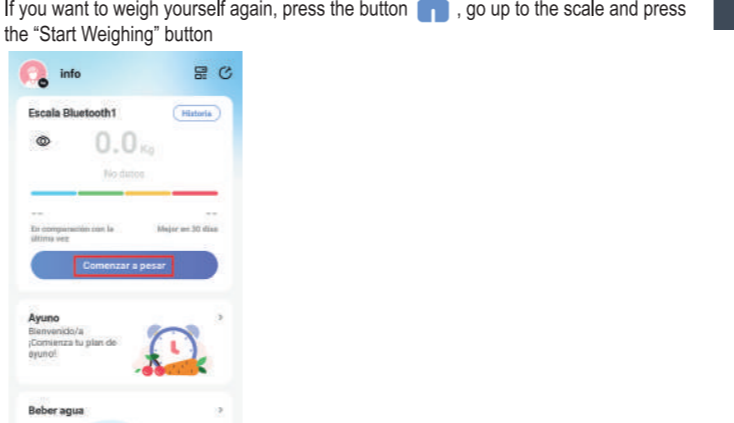
If you want to add new users to the scale press the button , then press "Add new user", complete the personal information and press the "Confirm" button
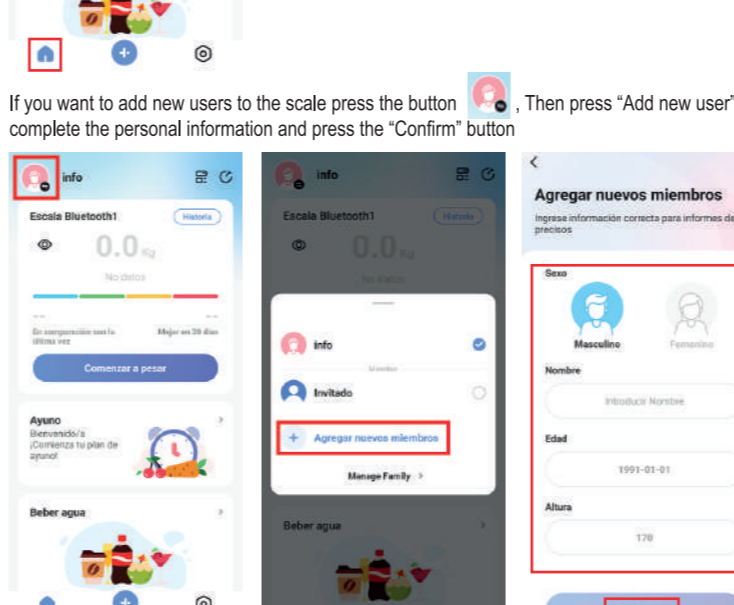
### Other functions:

The application incorporates new utilities such as:

- Fast**: Choose a meal plan and track your meals
- Drink water**: Monitor your water intake and set goals.
- Calorie**: Record the calories you eat and the weight of your food, as well as the goals you want to meet.
- Waist**: Record your waist and hip measurements
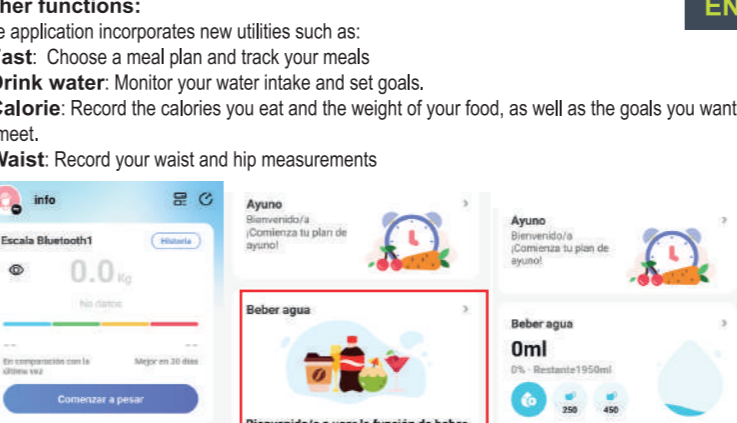
### Main Features

Capacity: 180 Kg / Power supply: 2X1.5V AAA battery / LCD display with backlight / 28x28 cm, 5 mm thickness tempered glass / Auto power on-off / Overload indication / Low battery indication / BT function connect to the mobile app / Measures weight, calories, body fat, muscle mass, bone mass and body water / APK-APP: oKOK International / BT frequency band: 2.4 GHz / Maximum emitted RF power < 100 mW

### Notes

- Keep scale surface dry to avoid accident by slipping.
- Please put the scale in a dry place to protect the internal electronic components.
- Please use micro cloth to clean the scale, not immersed in water.
- Please don't remove the scale except replace the battery.
- Please don't use heavy percussion or falls from a height.
- Please don't persistent put heavy on the scale to avoid damage to the high precision sensor.
- Please keep the scale away from wet floor, don't use the scale

- Please stand on the middle of scale in case of turnover.
- The platform of the scale is made of glass, please avoid sharp things scratching the surface, excessive forces or sudden shocks, unexpected fall-off will break the glass and injure the user.
- The Scale is not waterproof. Keep the surrounding environment dry, do not submerge it in water.
- use cloth with neutral cleaning liquids to clean the surface, keep away contacting with organic solvent nor chemicals.
- When the batteries are in low power status, please remove the batteries in case for leaking liquids could damage the scale.
- Please stand on the scale in the right direction.
- Please stand on the scale no less than 10 seconds if you measure body fat.
- Please confirm gender, height, weight data is correct.
- Low battery indication: 2.4 +/- 0.2V, display "Lo" and power off
- Overload indication: Overload more than 182.5 Kgs display "Err"

### Guarantee

**GUARANTEE**

DEVICE ITEM: \_\_\_\_\_  
Serial number: \_\_\_\_\_  
ESTABLISHMENT  
Purchase date: \_\_\_\_\_

**AQUATIS**  
+34 91 713 94 98  
INNOVA CELLULAR SL  
C/ Severo Ochoa N° 9 nave B8  
28521 Rivas-Vaciamadrid Madrid

**Note:** Do not expose the scale at extreme temperatures, the batteries could exploit. The scale is not waterproof. The BASCULA BT CON CALCULO DE LA COMPOSICION CORPORAL AQ1538 (BAS/AQ1538 (LCD26F)) has a 3-year warranty from the date of purchase. Defects caused by use and wear of the product do not come under warranty. Batteries are not covered by warranty. The scale only comes under warranty if it has not been tampered with and the defect it presents is due to factory reasons.

### Declaration of conformity

Through this INNOVA CELLULAR, it declares that the BASCULA BT CON CALCULO DE LA COMPOSICION CORPORAL AQ1538 (BAS/AQ1538 (LCD26F)) complies with the essential requirements and any other applicable or enforceable provisions of Directives 2014/30/EU, 2014/53/EU (RED), RoHS 2011/65/EU, annex II (EU) 2015/863. To see the complete declaration of conformity, you can access it through the following link: [https://www.innovacelular.com/descargas/declaracion/declaracion\\_basqa1538.pdf](https://www.innovacelular.com/descargas/declaracion/declaracion_basqa1538.pdf)

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### Safety warnings

**Safety and handling:** For more information on safety and handling, see the user manual. **Important Safety Information:** Handle the scale with care. They contain sensitive electronic components, including batteries, and may malfunction, affect performance, or cause damage if dropped, burned, punctured, crushed, disassembled, or exposed to excessive heat, liquid, or high-pressure environments, concentrations of industrial chemicals, including nearby evaporated liquefied gases such as helium. Do not use the scale if it is damaged. **Choking hazard:** The scale and the small parts included may present a choking hazard to small children or cause other injuries. Always keep them out of the reach of children.

### Other data

**Box contents:** Scale, 2 AAA batteries, user manual

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